

# DIGISPENSE® 1000\* CONTROLLED SUBSTANCE DISPENSING SYSTEM

## FEATURES AND BENEFITS

- High-speed brushless DC motor & controller
- Whisper-quiet motor
- Featuring Positive Displacement Ceramic Pumps
- Easy to control from the USB or RS232 interface
- Factory-made tubing kits with bottle cap connections
- Precision volumetric flask for accuracy checks
- Factory pre-calibrated motor/base pump module
- Optional all-in-one Controlled Substance Dispensing Station: portable, secure & easy to clean\*\*
- CE certified

The IVEK DigiSpense® 1000\* is a precise, high speed, low volume, positive displacement liquid dispensing system. The DigiSpense® 1000 is capable of extremely high accuracy and repeatability and is often used for dosing pharmaceuticals, where inventory control and accuracy are critical. The wetted components include 99.8% pure Alumina ceramic piston/cylinder modules. These components are mated together to a clearance of 120 millionths of an inch. With only one moving part in the fluid path, systems operate for many years without measurable wear. A USB and RS232 interface port on the controller offers our customers a turnkey solution supported by most addiction treatment software.

IVEK Corporation has been manufacturing precision positive displacement dispensing systems specifically for use in addiction management programs for over 30 years. These systems are noted for their reliability, low maintenance and high precision fluid handling and are placed in treatment clinics worldwide.

\*Enhanced replacement for IVEK's DIGISPENSE 10/RS232


\*\*See page 2 for more information.



DigiSpense® 1000 Controller\*  
28.3cm x 14.3cm x 45.2cm



single-end motor base module  
7cm x 12.8cm x 10cm

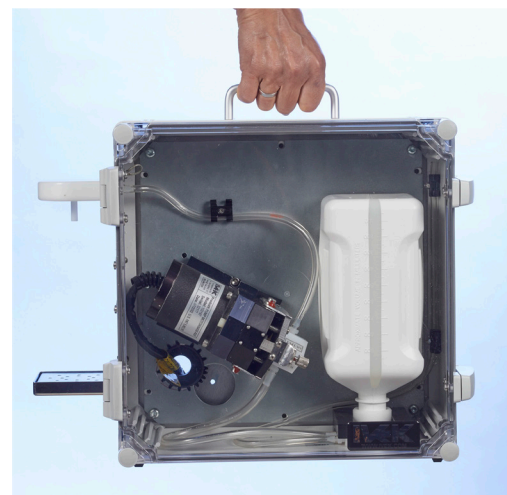

all-in-one portable, secure, easy to clean,  
controlled substance dispensing station option available\*\*  
14cm x 39.5cm x 35cm





10 FAIRBANKS ROAD NORTH SPRINGFIELD VERMONT USA 05150  
TEL (01) 802-886-2238 FAX (01) 802-886-8274  
TOLL FREE IN NORTH AMERICA 800-356-4746




**CSDS\*****CONTROLLED SUBSTANCE DISPENSING STATION\***

**OPTIONAL ALL-IN-ONE  
CONTROLLED SUBSTANCE  
DISPENSING STATION:  
PORTABLE, SECURE & EASY TO CLEAN!**

- Featuring IVEK's Positive Displacement Ceramic Pumps
- Precision Volumetric Flask for Accuracy Checks
- Factory Pre-Calibrated Motor/Base Pump Module
- Shorter Tubing Runs for Faster Priming & Less Waste
- Fixed Dispense Tip for No Mess Dosing
- Adjustable Spill Tray Locations for Varying Cup Sizes
- Lockable Latches (locks not supplied)
- Sealed Front Cover allowing for up to 1L of Fluid Spill inside Enclosure without Leaks
- Zero-Transfer Bottle Holder & is Compatible with All Major Medicine Manufacturers
- Compatible with Existing IVEK DigiSpense 10/DC/1000 Pumps
- FREE Installation & Calibration of Existing Pump with CSDS Purchase!
- CSDS with Pump may be returned to IVEK for Standard Cleaning & Calibration

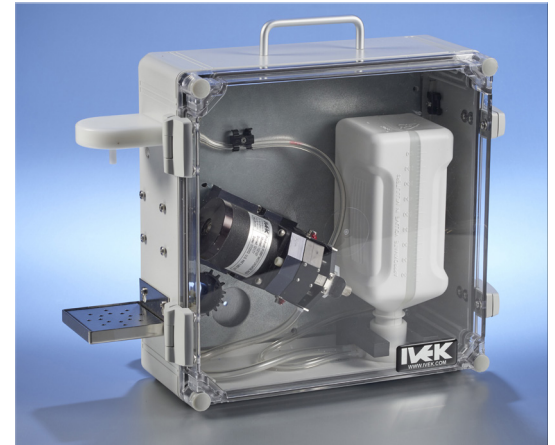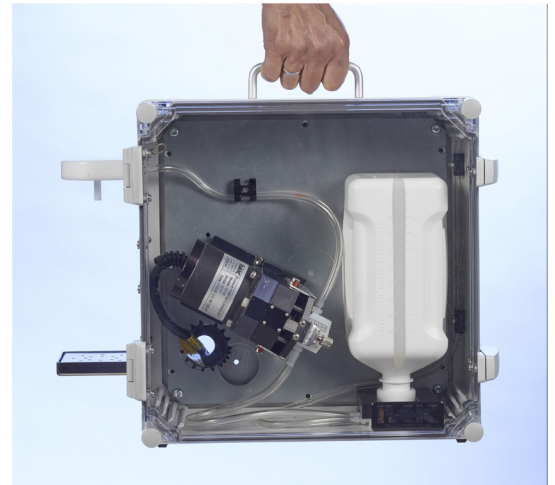

14cm x 39.5cm x 35cm

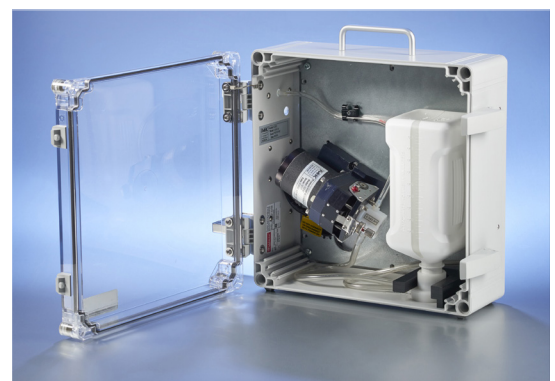

\*Optional With IVEK's DigiSpense® 10/RS232 & DigiSpense® 1000

**PRECISION ENGINEERED FLUIDIC SOLUTIONS SINCE 1979**

IVEK Corporation is committed to quality and customer support. Our product has always been relied upon for critical applications and our policy states customer satisfaction as the number one priority. As our customer you will experience exceptional performance and reliability for your most demanding applications. IVEK Corporation is ISO 9001 and ISO 13485 certified.



IVEK CORPORATION • 10 FAIRBANKS ROAD • NORTH SPRINGFIELD • VERMONT • 05150 • USA  
TEL: (01) 802-886-2238 • FAX: (01) 802-886-8274 • TOLL FREE IN NORTH AMERICA: 800-356-4746